

## Contents

<b>Checklist Protocol 1</b> .....	1
Sensitivity documentation table (Rectangular vape) .....	2
Suggested Friends Device Placement: .....	3
<b>Checklist Protocol 2</b> .....	4
VM-FRIENDS Test Result .....	4

## Checklist Protocol 1


Name: Md Rafi Islam		Vape Brand: Juul		Date: 04/25/2024	
Task no.	Steps			Status (Y/N?)	
1	Charge the vape.			Y	
2	Charge the FRIENDS device.			Y	
3	Mark RF sensor location on the FRIENDS enclosure				
4	Draw a 7mm grid on the vape			Y	
5	Wrap the vape's mouthpiece with putty.				
6	Attach the vape to the vacuum line.			Y	
7	Identifying the sensitive zones	7a	Select a side (front, back, right, left, or bottom), place the FRIENDS on top left corner of that side.	Y	
		7b	Initiate puff with the vacuum and very slowly slide the FRIENDS from top to the bottom. <b>Make sure the vape is puffing by observing the puff indicator light on the vape.</b>	Y	
		7c	Observe the Green LED of the FRIENDS	For stable green, mark the location with "X" in the grid on the vape.	Y
				For no LED illumination, left the space blank in the grid.	Y
		7d	After marking the sensitive spot on the vape, start sliding the FRIENDS device again and complete the whole row.		Y
		7e	Then place the FRIENDS on the middle of the vape and repeat steps 7b to 7d.		Y
		7f	After finishing scanning this row with FRIENDS, go to the next row and so on.		Y
		7g	In this way, the front side of the vape will be scanned for sensitivity zones and the grid on the vape will be marked with the "X" symbol if the vape is sensitive to FRIENDS device.		Y
		7h	Then place the FRIENDS on another side (e.g. front, back, right, left, or bottom) of the vape and start scanning again. Repeat steps 7a to 7e and populate the documentation table.		Y
		7i	Finally, after completing scanning all the sides of the vape, take a picture of each side of the vape and insert the images in the checklist document table below.		Y
		7j	If there are multiple sensitive zones, select a possible FRIENDS mounting location based on signal strength and usability.		Y
7k	Take a picture and insert it at the end of the checklist documentation.		Y		

## Sensitivity documentation table (Rectangular vape)

For a stable green LED, mark the location with the letter "S" in the documentation



Front view

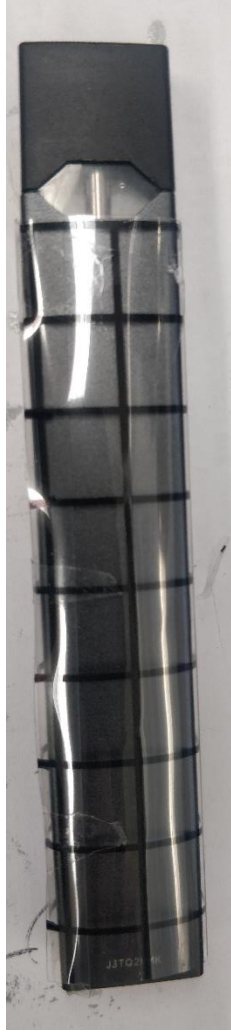

Back view




Left view



Right view

Suggested Friends Device Placement:

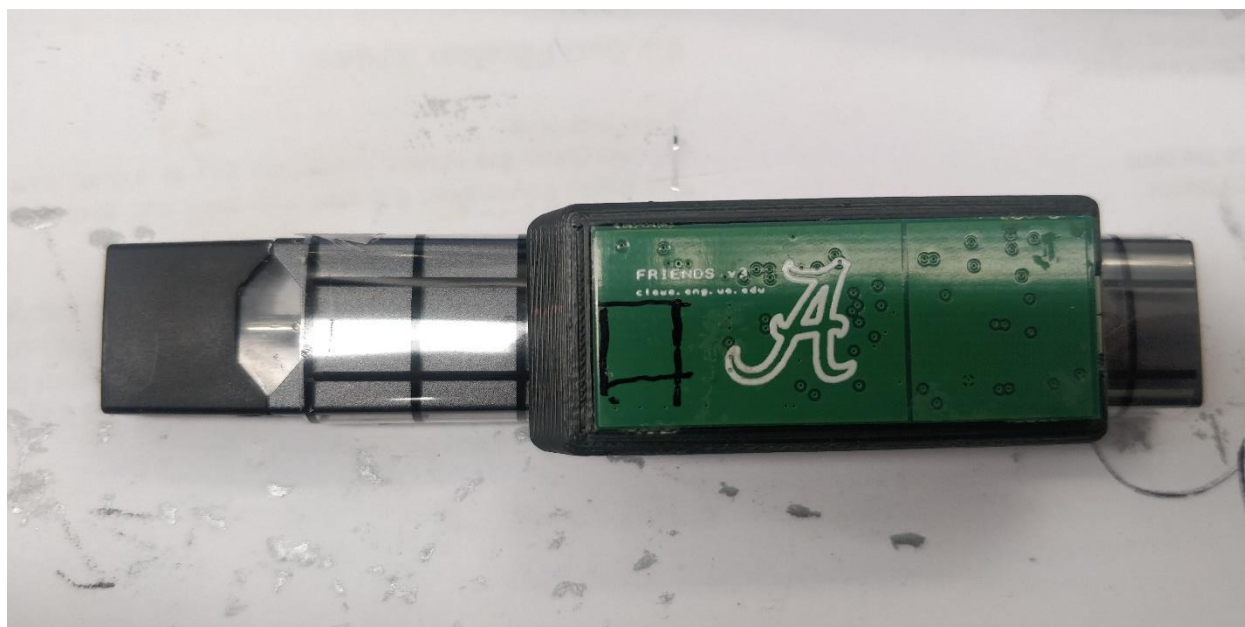

RF sensor Location: Front center

Note: The side with the JUUL writing is the front side of the vape

## Checklist Protocol 2

Task no.	Steps		Status (Y/N?)	
1	Charge the vape.		Y	
2	Charge the FRIENDS device.		Y	
3	Wrap the vape's mouthpiece with putty.		Y	
4	Attach the vape to the vacuum line.		Y	
5	Vaping Machine for laboratory testing of FRIENDS sensors (VM-FRIENDS) Testing	Y	Recharge the FRIENDS device and the vape.	Y
		Y	Attach the vape with a vacuum line leading to the Vaping Machine for laboratory testing of FRIENDS sensors (VM-FRIENDS)	Y
		Y	Open VM-FRIENDS GUI and Set the following parameter:  Select input parameters: Volume and Duration Volume (mL): 120 Inhale Flowrate (mL/s): 40.0 Exhale Flowrate (mL/s): 60.0 Inhale Time (s): 3 Exhale Time (s): 2 Interpuff Time (s): 18 Number of Puffs: 80	Y
		Y	Click on Start Puffing	Y
		Y	Read and store the FRIENDS device data.	Y
		Y	After 80 puffs, repeat steps 8a-8e till the vape runs out of vape juice.	Y
		Y	After completing step 8f, fill out the table given at the end of the document.	Y

### VM-FRIENDS Test Result

Vape's Total Puff Capacity	FRIENDS Total Puffs Registered	Thermistor Used for Testing	Thermistor Needed for Operation
185	185	No	No